


Part designed for 3D Printed Prototype Evaluation
Original Design and Testing Procedure Found in
Army Report TM-68-30

Port A: FIO2 Supply

Port B: To Patient, Breathing valve is required

Port C: Exhaust, Fit with N95 Filter, can be any
throttling valve for evaluation

B

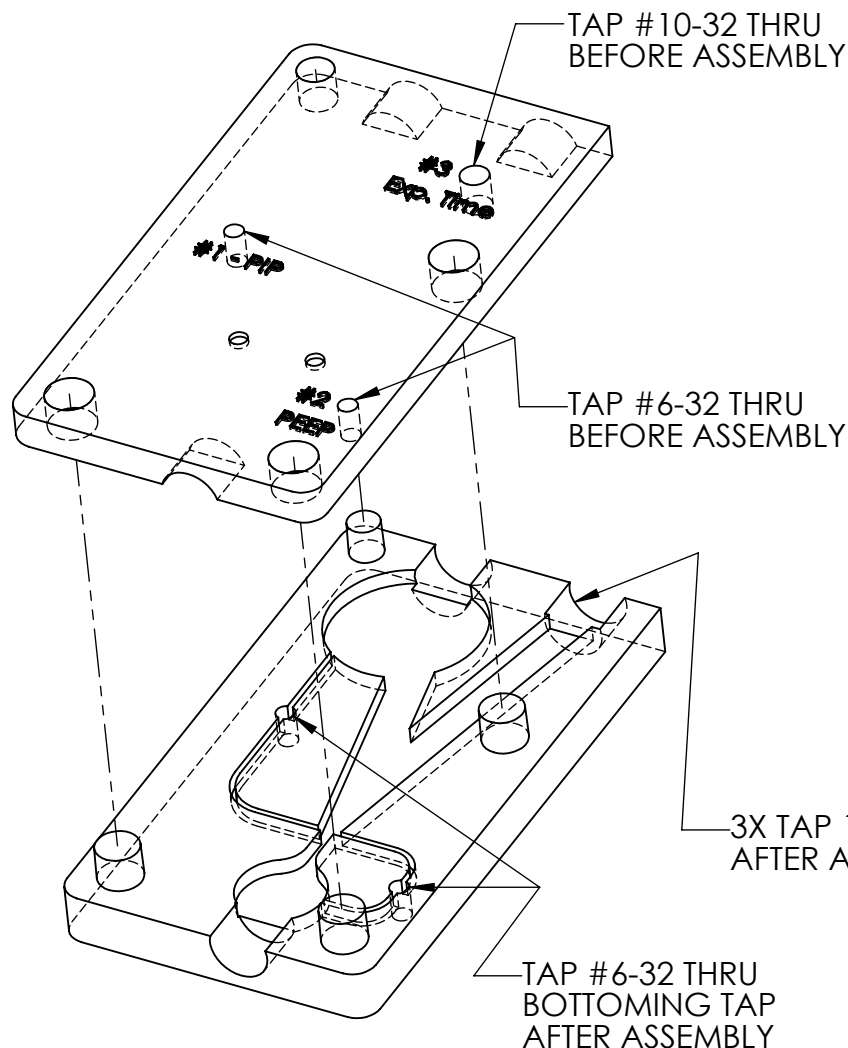

A


Indents mark location of
control channels. These can
be drilled in lid before
assembly to install pressure
sensors. Not currently needed

Calibration Procedure (See Army
report for background):

1. Set input pressure
2. Adjust Screw #1 to set PIP, the
pressure at which the device
switches to expiration
3. Adjust Screw #2 to set the
pressure at which the device
switches back to inhalation.
This is the PEEP.
4. Adjust Screw #3 to set the
duration of the expiration
phase.

It is assumed that gauges can be
attached to ports as required.

Target values for PEEP, PIP, Tidal
Volume, Respiratory Rate are
documented elsewhere

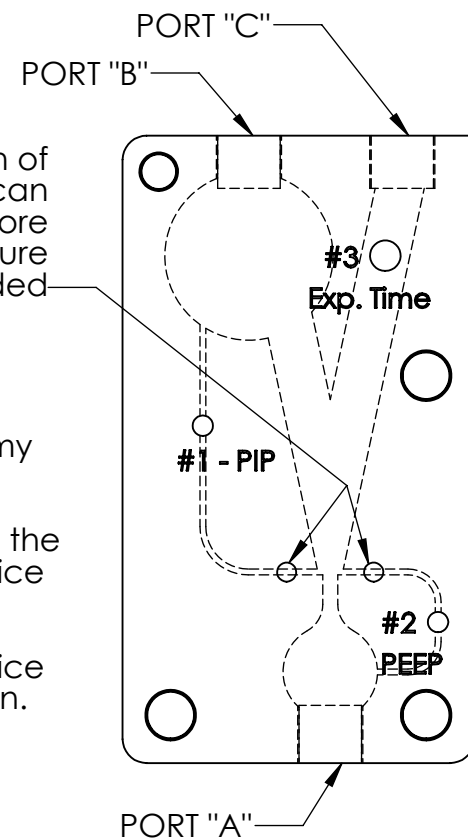

B

A

Ryan Whitney
m.ryan.whitney@gmail.com

Army Ventilator w/PEEP For 3D Print Prototype Rev 02

SIZE	DWG. NO.	REV
A	Army Vent PEEP 02, Assm	
SCALE: 1:1	WEIGHT: 0.101	SHEET 1 OF 1